

Aidite[®]

A-IS Pro 3 Digital Scanner



A-IS Pro 3

Aidite A-IS Pro 3 is a 3D dental scanner with many powerful functions for multiple applications covering scans of impressions, plaster models, articulators, implant abutments, etc. Its high-resolution cameras, advanced algorithms and user-friendly interface provide you with superior quality data, providing optimal results.



Aidite

Aidite A-IS Pro3

Advanced Desktop Scanner




Ultra-high Accuracy

Scan accuracy less than $8\mu\text{m}$ assures reliable scan result, satisfying the requirement of various applications in the dental industry.



Unparalleled Details

Equipped with dual 5.0 MP high resolution cameras, A-IS Pro 3 can capture every detail of the components meticulously.

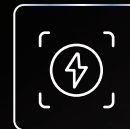

Fast Speed

Bite scan: 7s
Upper/Lower jaw scan: 14s
1~4 dies scan: 18s
5~8 dies scan: 30s
Impression scan: 66s

Powerful Software

Flexible Workflow


Set up your own workflow as you prefer, increasing work efficiency like never before, all-in-one and triple tray are available to get complete scan data with less steps.

Backstage Data Processing

Scan and backstage data processing can be carried out simultaneously, making the whole process more efficient.

Flexible Workflow


The new software associated with A-IS Pro 3 can support versatile scan modes, including Multi Path Scan, High Quality Scan, AI Scan, Texture Scan, HDR Mode Scan, Reduce High Brightness Scan, etc.

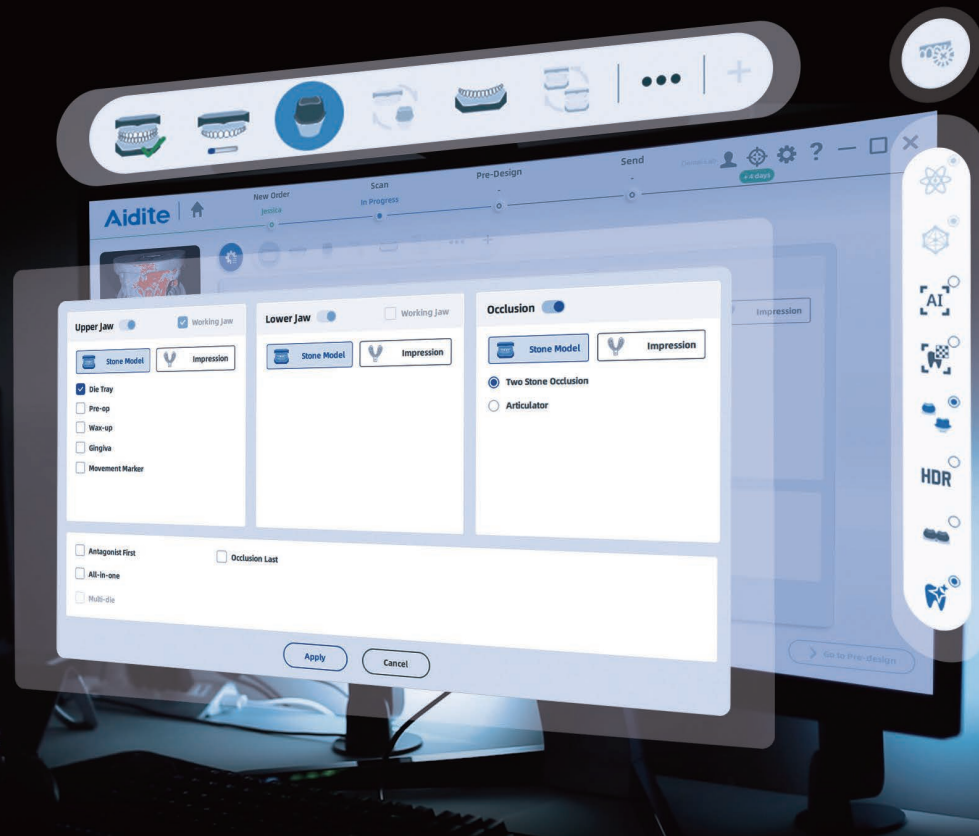

Applications



Implant Abutment

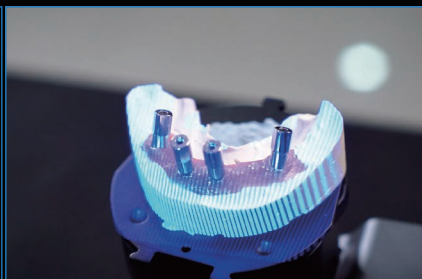

Implant Bridge



Scanbody




Texture



Model



Articulator

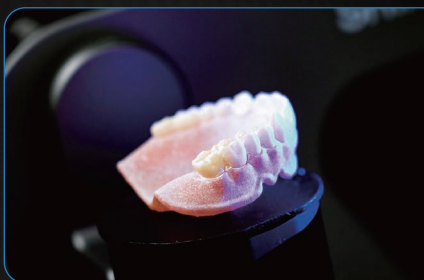

Full Denture



Impression



Post and Core

Technical Specifications

Model	A-IS Pro 3
Camera Resolution	2*5.0 MP
Accuracy	≤ 8 μm
Scan Time	Bite scan: 7s; Upper/Lower jaw scan: 14s; 1~4 dies scan: 18s; 5~8 dies scan: 30s; Impression scan: 66s
Articulator Scan	Support articulator transfer
Output Format	STL, PLY, OBJ
Light Source	Blue Light
Dimension	385*260*270 mm
Weight	5.0 Kg
Interface	USB 3.0
Power Supply	DC 24V

Aidite[®]

Aidite (Qinhuangdao) Technology Co.,Ltd.

Tel : 0086-335-8587898

Web: www.aidite.com

Email: info@aidite.com



www.aidite.com